

EXTERNAL LUBRICATION PACK

Spray Magnesium Stearate on your STYL'ONE

STYL'One




ME-050_EN_C (07 2019)

EXTERNAL LUBRICATION PACK

STYL'One



5g capacity container

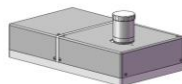

Micro-dosing unit



Pulsed air blow cabinet

- The micro-dosing unit has been specifically designed for single punch presses. Container can hold up to 5 g of Magnesium Stearate.
- A puff of Magnesium Stearate is only generated when the spraying nozzle is above the die.
- Upper and lower punches are lubricated at the same time.
- The press remains clean since only the minimum amount of lubricant is sprayed onto the tooling
- Adjustable pulsed air blow from 0 to 5 bars



Micro-dosing unit + pulsed air blow cabinet connected to a STYL'ONE Tableting Instrument

Technical data	Micro-dosing Unit	Pulsed Air blow cabinet
Inlet compressed air requirement	-	6 bars
Adjustable pulsed air blow	-	0 to 6 bars
Lubricant container Capacity*	5 g	-
Standard voltage	-	Powered by STYL'One (24 VDC)
Weight (kg / lbs.)	3.5 kg (7.7 lbs.)	21 kg (46 lbs.)
Dimension (length / width / height)	381mm x 221mm x 590mm 15" x 8" 7/10 x 23" 2/10	
Swivels wheels	-	4 (2 with brakes)

* Capacity is product dependant. Indivative value based on Peter Greven LIGAMED MF-3-V grade of vegetable Magnesium Sterate

